



# microTargeting™ Solution for Monteris® Laser Ablation



**Patient-specific stereotactic system that allows a surgeon to quickly and accurately place a Monteris® mini-bolt used to support the Monteris® Robotic Probe Driver**

The microTargeting™ solution for Monteris® laser ablation is one of the latest technology advances introduced as part of FHC's Connected Care product portfolio and our commitment to working with the medical community in improving patient care. The adapter is designed to act as the interface between the microTargeting Platform and the Monteris® NeuroBlate® system for MRI-guided neurosurgical ablations and lesions.

The adapter is attached to a standard microTargeting Platform via a mounting hub and secured with provided thumbknobs. Once the adapter is correctly oriented in the hub, the Monteris® NeuroBlate® Guide Rail and laser positioner may be attached. Both entry and target mounting hubs attach to the adapter and allow a change in planned trajectory.

microTargeting™ is a trademark of FHC, Inc.

- ✓ **Flexibility:** allows the Monteris® NeuroBlate® system to be used with FHC's patient-customized frameless stereotactic platform for neurosurgical procedures
- ✓ **MR Safe:** allows use in MRI environment as well as the OR
- ✓ **Adaptable:** attaches directly to the Monteris® NeuroBlate® system and to the microTargeting Platform
- ✓ **Steam Sterilizable:** for rapid turn around and ease of use
- ✓ **Convenience:** can be sterilized with the WayPoint™ Sterilization Kit in the SteriSuite Case
- ✓ **Unique:** target offset and entry offset mounting hubs allow change of the planned trajectory for a given microTargeting™ Platform
- ✓ **Integrated:** for epilepsy cases, reuse anchors and scans used to implant SEEG depth electrodes

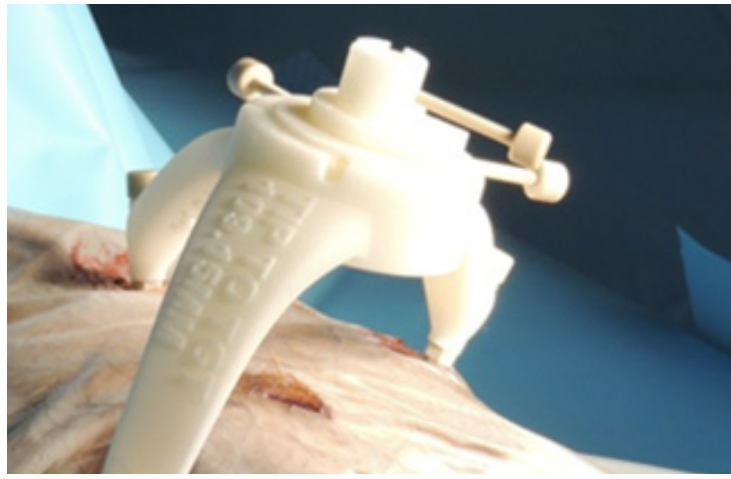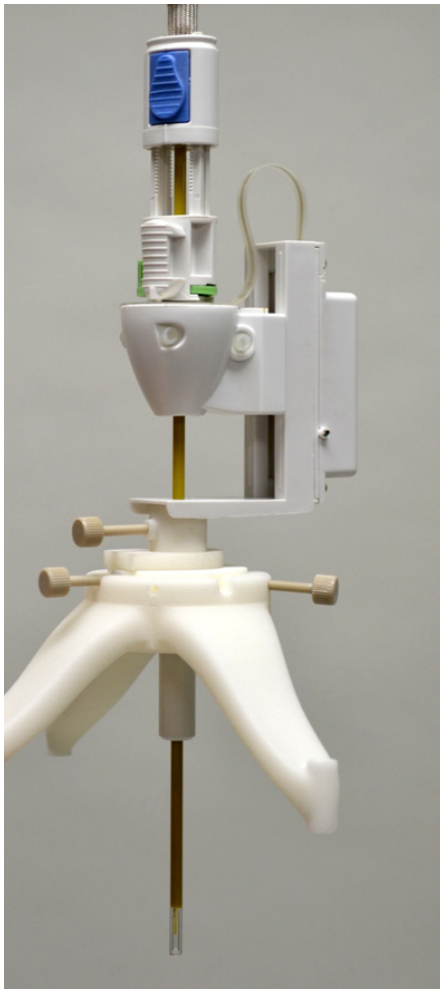

### Ordering Information:

66-WP-NB – Monteris® NeuroBlate®  
microTargeting Platform Adapter

### Designed for use with:

NBD-01 – NeuroBlade Advanced Probe Driver

### Additional Items Required for Operation:

MP-KIT-XX: microTargeting Platform\*

\* Platform includes disposable implant kit and disposable surgical kit

L022-03 (Rev B0 2018-02-01)



#### FHC, Inc.

1201 Main Street  
Bowdoin, ME 04287 USA  
Fax: +1-207-666-8292  
www.fh-co.com



#### FHC Europe

(TERMOBIT PROD srl)  
42A Barbu Vacarescu Str, 3rd Fl  
Bucharest 020281 Sector 2  
Romania

#### FHC Latin America

Calle 6 Sur Cra 43 A-200  
Edificio LUGO Oficina 1406  
Medellín-Colombia